


JUPITER

PEEK Cage for Posterior Lumbar Interbody Fusion (PLIF)

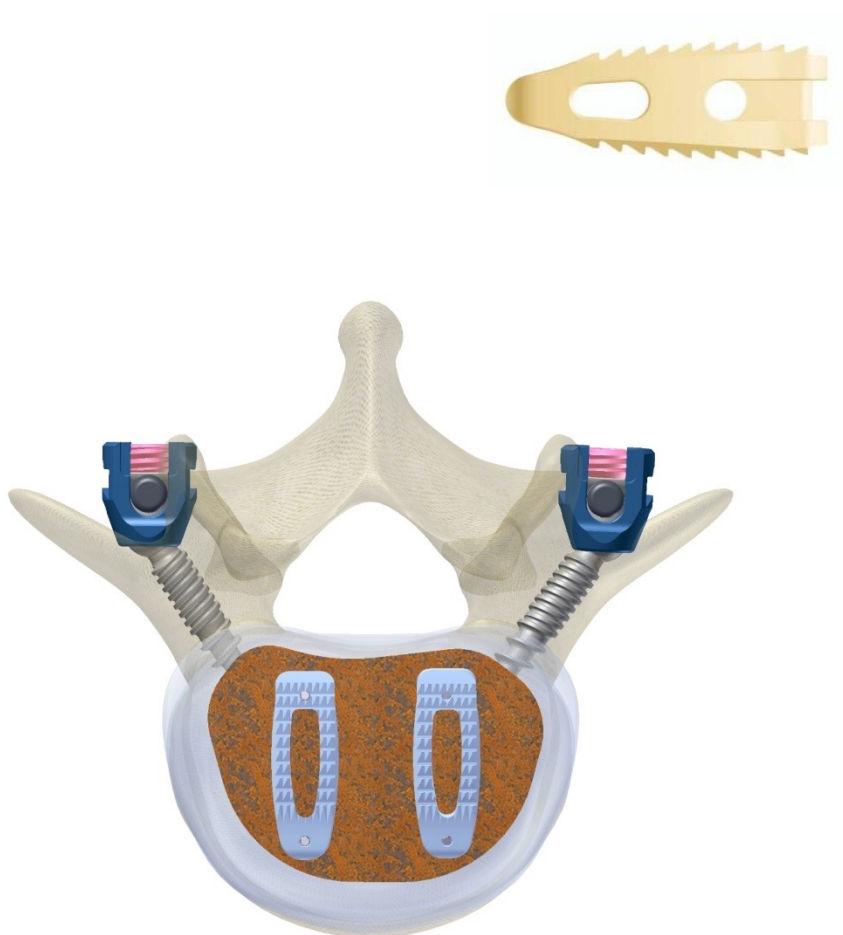

SYSTEM OVERVIEW

The elegant simplicity of the JUPITER™ design along with the flexibility of its use in a posterior approach appeals to surgeons that prefer a consistent and reliable implant. The JUPITER™ PLIF PEEK Cage System was specifically designed for surgeons that prefer a traditional posterior approach for bilateral implant placement. The specific shape and the pure PEEK composition of this implant ensure excellent bone fusion.

Implant Features — PEEK-OPTIMA®

Safety

- Elasticity prevents sinking: less risk of endplate penetration
- Withstands repeated autoclave sterilizations

Medical Imaging (X-ray, CT, MRI)

- Radiolucent material to allow accurate fusion follow-up
- Embedded gold markers to facilitate implant placement verification

Biocompatibility


- No carbon fibres: no risk of inflammatory reaction
- CE-marked and cleared by FDA for long-term implantation in the human body

Enhanced Bone Fusion

- No stress shielding
- Optimum load repartition
- Bone growth enhancement thanks to micro-movements

Modulus of Elasticity

- Elastic Modulus between cancellous and cortical bone
- Ideal load sharing implant

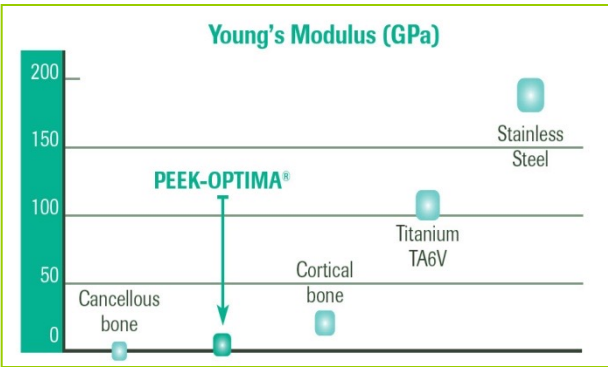

Implant Features — Cage Design

- Bulleted nose minimizes insertion forces and provides an easy, atraumatic and self-distracting insertion.
- Convex-shaped implants are designed to fit patient anatomy and to allow more accurate sizing.
- X-ray markers facilitate radiologic identification.
- Biocompatible PEEK Optima allows the growth of the bone in the cage to be visualized.
- Sharp teeth on the surface of the implant ensure primary stability and prevent migration of the cage.



Indications

The intended purpose of the JUPITER™ PLIF PEEK Cage is that it be used for interbody fusion and/or resection or excision of the vertebral bodies of the thoracolumbar spine.

JUPITER™ PLIF PEEK Cage is indicated for the treatment of degenerative disc disease, disc herniation, foraminal stenosis, tumor, trauma, deformity (including scoliosis, spondylolisthesis and retrolisthesis) and failed previous fusion. JUPITER™ was designed to be placed through a posterior incision.

NOTE:

Since the JUPITER™ PLIF PEEK Cages were not developed as “stand-alone” implants, the use of additional posterior instrumentation (for example with NEPTUNE™ Pedicle Screws) is strongly advised.


IMPLANT INFORMATION

JUPITER™ PLIF PEEK Cage, Footprint 10 X 28 mm

HEIGHT (MM)	REF
8	MOI 47005008
9	MOI 47005009
10	MOI 47005010
11	MOI 47005011
12	MOI 47005012
13	MOI 47005013
15	MOI 47005015
17	MOI 47005017

JUPITER™ PLIF PEEK Cage, Footprint 10 X 24 mm

HEIGHT (MM)	REF
8	MOI 47030008
9	MOI 47030009
10	MOI 47030010
11	MOI 47030011
12	MOI 47030012
13	MOI 47030013
15	MOI 47030015
17	MOI 47030017

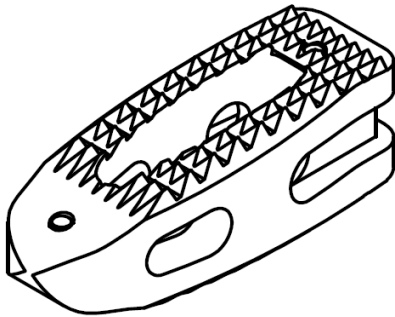



USA OFFICE

2404 Hollywood Blvd.
Hollywood, FL 33020

PUERTO RICO OFFICE

1660 Santa Ana Street
Esq. Fidalgo Diaz,
San Juan, Puerto Rico 00909

E-mail: info@madisonorthoinc.com
Website: www.madisonorthoinc.com

