

PRIMA™

Thoracolumbar Spine Locking Plate System

System Overview

The PRIMA™ Thoracolumbar Spine Locking Plate (TSLP) is a versatile, low-profile plating system indicated for anterior stabilization of the thoracic and lumbar spine (T1–L5). The system provides a 4.5 mm construct profile, plate lengths from 40 mm to 109 mm in 3 mm increments, 5.5 mm cancellous locking screws, and instrumentation designed to simplify construct placement.


Implant Features

■ TSLC Plates

- Comprehensive line of plates for mono- and bisegmental instrumentations
- Low profile of 4.5 mm with a width of 26 mm
- Range of plate lengths (in 3-mm gradations)
- Anatomical plate curvature
- Screws lock with angular stability
- Alternative screw hole selection
- Central holes for bone graft visualisation
- Titanium alloy (Ti-6Al-4V)


■ TSLC Screws

- Dual-core design
- Self-tapping
- Screw-head thread provides a form-fit with the plate
- Proximal cortical thread and distal spongiosa thread for appropriate anchoring in the bone
- A constant outer diameter increases resistance to tearing out
- Blunt screw tip for greater safety

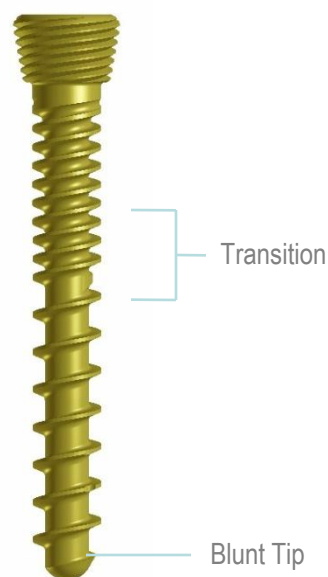

Indications for Use

The PRIMA™ TSLP System is indicated for use via the lateral or anterolateral surgical approach in the treatment of thoracic and thoracolumbar (T1– L5) spine instability resulting from fracture (including dislocation and subluxation), tumor, degenerative disc disease (defined as back pain of discogenic origin with degeneration of the disc confirmed by patient history and radiographic studies), scoliosis, kyphosis, lordosis, spinal stenosis, or a failed previous spine surgery.

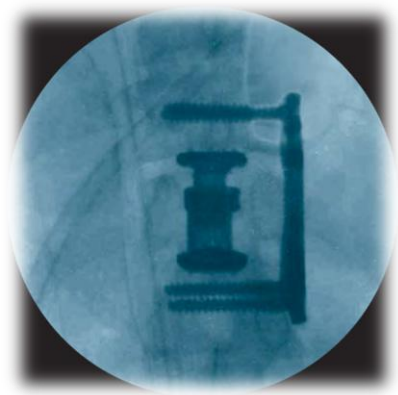

Implant Information [Material in Titanium Alloy (Ti-6Al-4V)]



TSLP Plate, 5.5

REF (TIT)	LENGTH
MOI 37020040	40 mm
MOI 37020043	43 mm
MOI 37020046	46 mm
MOI 37020049	49 mm
MOI 37020052	52 mm
MOI 37020055	55 mm
MOI 37020058	58 mm
MOI 37020061	61 mm
MOI 37020064	64 mm
MOI 37020067	67 mm
MOI 37020070	70 mm
MOI 37020073	73 mm
MOI 37020076	76 mm
MOI 37020079	79 mm
MOI 37020082	82 mm
MOI 37020085	85 mm
MOI 37020088	88 mm
MOI 37020091	91 mm
MOI 37020094	94 mm
MOI 37020097	97 mm
MOI 37020100	100 mm
MOI 37020103	103 mm
MOI 37020106	106 mm
MOI 37020109	109 mm



TSLP Screws, 5.5


REF (TIT)	LENGTH
MOI 37021022	22 mm
MOI 37021024	24 mm
MOI 37021026	26 mm
MOI 37021028	28 mm
MOI 37021030	30 mm
MOI 37021032	32 mm
MOI 37021034	34 mm
MOI 37021036	36 mm
MOI 37021038	38 mm
MOI 37021040	40 mm
MOI 37021042	42 mm
MOI 37021044	44 mm
MOI 37021046	46 mm
MOI 37021048	48 mm
MOI 37021050	50 mm
MOI 37021052	52 mm
MOI 37021054	54 mm




Bone Graft Fixation Screws, 3.0 (Optional)

REF (TIT)	LENGTH
MOI 37022030	30 mm
MOI 37022032	32 mm
MOI 37022034	34 mm
MOI 37022036	36 mm
MOI 37022038	38 mm
MOI 37022040	40 mm
MOI 37022042	42 mm
MOI 37022044	44 mm
MOI 37022046	46 mm
MOI 37022048	48 mm
MOI 37022050	50 mm

Instrument Set Details

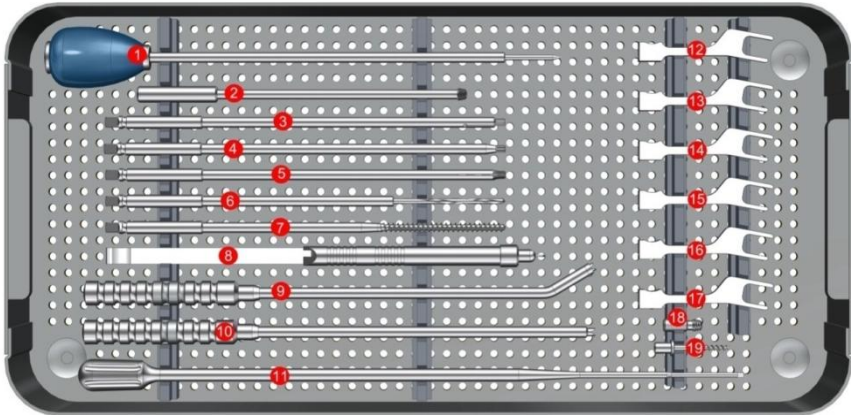

REF CODE	DESCRIPTION	UNITS
1 3100016	Awl, Ø3.3 mm	1
2 3100015	TSLP Plate Holder	1
3 3100023	Tightening Screwdriver Shaft SW4.5	1
4 3100022	Selfholding Screwdriver Shaft SW4.5	1
5 3100021	Driver SW4.0, for Threaded Drill Guide & Fixation Pin	1
6 3100018	Drill Bit, Ø3.3mm	2
7 3100017	Tap for Cancellous Bone Screws Ø5.5mm	1
8 1700008	Depth Gauge	1
9 3100014	Forked Impactor, Angled	1
10 3100013	Forked Impactor, Straight	1
11 3100010	Screwdriver SW2.5, for VBR Implant	1
12 3100001	Tips for Parallel Distractor, 12L	1
13 3100002	Tips for Parallel Distractor, 12R	1
14 3100003	Tips for Parallel Distractor, 15L	1
15 3100004	Tips for Parallel Distractor, 15R	1
16 3100005	Tips for Parallel Distractor, 20L	1
17 3100006	Tips for Parallel Distractor, 20R	1
18 3100020	Threaded Drill Guide	4
19 3100019	Temporary Fixation Pin	2

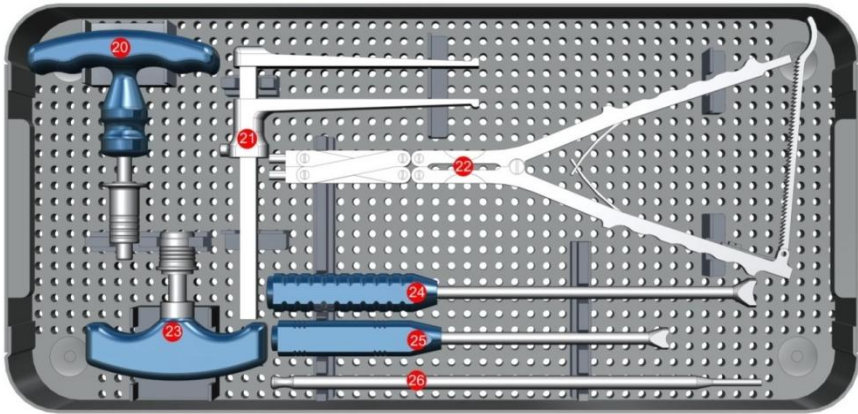

REF CODE	DESCRIPTION	UNITS
20 3100025	Torque-limiting T-Handle, 7NM	1
21 3100007	Inter Vertebral Measuring Device	1
22 3100009	Parallel Distractor for VBR Implant	1
23 3100024	Quick Coupling T-Handle	1
24 3100011	Impactor, Crescent, Large	1
25 3100012	Impactor, Crescent, Small	1
26 3100008	Graft Pusher	1