


AMAZON  **LCP** TM

LCP Hook Plate 3.5 mm

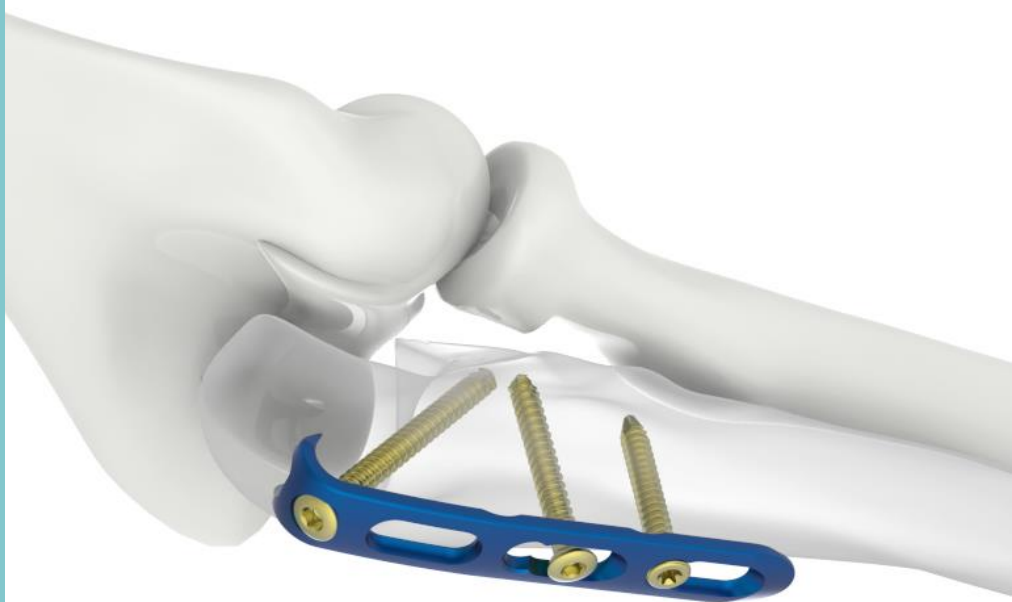


INTRODUCTION

The AMAZON™ LCP Hook Plate 3.5 mm is a locking compression plate that merges locking screw technology with conventional plating techniques. The plate features pointed hooks to aid in plate placement and holes for both lag screws and locking screws. Elongated LCP holes make the plate and its fixation more flexible and allow for controlled compression. The spring effect facilitates reduction and a stable tension-band technique.

The AMAZON™ LCP Hook Plate 3.5 mm enables tension-band plating with use of 3.5 mm screws. The hook plate was primarily intended for the treatment of simple fractures of the olecranon as well as osteotomies of the olecranon following distal humeral fracture. It has also shown benefits in avulsion fracture of the distal tibia and fibula. A single 3-hole plate can be used in multiple locations, applied to either the right or left side of the anatomy, resulting in less inventory required.

Implant Features

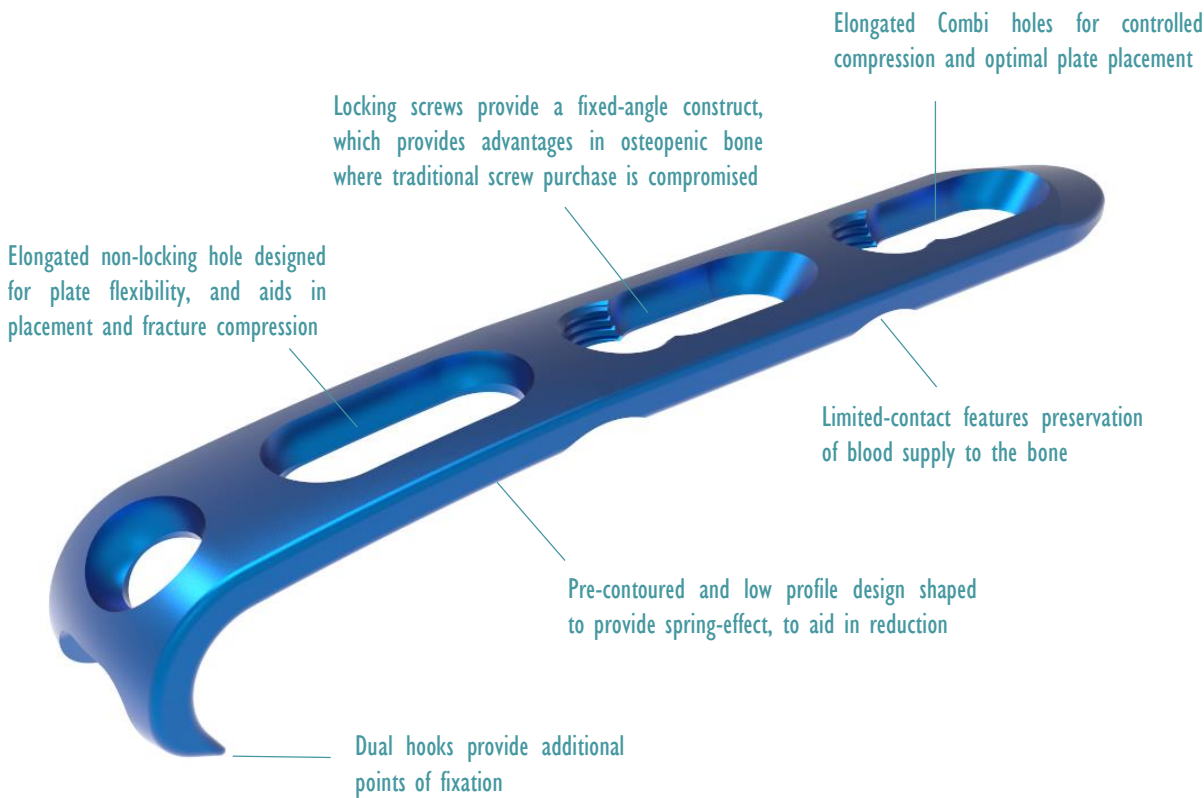

Indications

The AMAZON™ LCP Hook Plate 3.5 mm is indicated for fractures, osteotomies and non-unions of small bones, including the ulna, radius, tibia and fibula, particularly in osteopenic bone.



IMPLANT INFORMATION




LCP Hook Plate 3.5 mm

REF (TIT)	HOLES
MOI 35296003	3



3.5 mm Locking Screws

REF (TIT)	LENGTH
MOI 35228012	12 mm
MOI 35228014	14 mm
MOI 35228016	16 mm
MOI 35228018	18 mm
MOI 35228020	20 mm
MOI 35228022	22 mm
MOI 35228024	24 mm
MOI 35228026	26 mm
MOI 35228028	28 mm
MOI 35228030	30 mm
MOI 35228032	32 mm
MOI 35228034	34 mm
MOI 35228036	36 mm
MOI 35228038	38 mm
MOI 35228040	40 mm
MOI 35228042	42 mm
MOI 35228044	44 mm
MOI 35228046	46 mm
MOI 35228048	48 mm
MOI 35228050	50 mm
MOI 35228052	52 mm
MOI 35228054	54 mm
MOI 35228056	56 mm
MOI 35228058	58 mm
MOI 35228060	60 mm

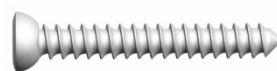

3.5 mm Cortex Screws

REF (TIT)	LENGTH
MOI 33500014	14 mm
MOI 33500016	16 mm
MOI 33500018	18 mm
MOI 33500020	20 mm
MOI 33500022	22 mm
MOI 33500024	24 mm
MOI 33500026	26 mm
MOI 33500028	28 mm
MOI 33500030	30 mm
MOI 33500032	32 mm
MOI 33500034	34 mm
MOI 33500036	36 mm
MOI 33500038	38 mm
MOI 33500040	40 mm
MOI 33500042	42 mm
MOI 33500044	44 mm
MOI 33500046	46 mm
MOI 33500048	48 mm
MOI 33500050	50 mm
MOI 33500052	52 mm
MOI 33500054	54 mm
MOI 33500056	56 mm
MOI 33500058	58 mm
MOI 33500060	60 mm



| TRAUMA

USA OFFICE

2711 S Ocean Dr.
Unit 2704
Hollywood, FL 33019

PUERTO RICO OFFICE

1660 Santa Ana Street
Esq. Fidalgo Diaz,
San Juan, Puerto Rico 00909

E-mail: info@madisonorthoinc.com
Website: www.madisonorthoinc.com

